FINEVISION TRIFOCAL OPTIC



FINEVISION

The first trifocal diffractive IOL

PhysIOL

Technical specifications

Commercial name	Micro F		Pod F	
Material	25% hydrophilic acrylic		26% hydrophilic acrylic	
Overall diameter	10.75 mm		11.40 mm	
Optic diameter	6.15 mm		6.00 mm	
Optic	Biconvex aspheric (-0.11µSA) trifocal diffractive FineVision			
Filtration	UV & blue light			
Refractive index	1.46			
Abbenumber	58			
Angulation	5°			
Additional power	+ 1.75D for intermediate vision & + 3.50D for near vision			
Injection system	Medicel Viscoject Bio 1.8 / 2.2 Medicel Accuject 1.8 / 2.0 /2.1 / 2.2		Medicel Accuject 2.0 from 6D to 24.5D Medicel Accuject 2.1 / 2.2 from 25D to 35D	
Incision size	≥ 1.8 mm		≥ 2.0 mm	
Spherical power	10D to 35D (0.5D steps)		6D to 35D (0.5D steps)	
Square edge		36	;0°	
Nominal manufacturer A constant	118.80		118.95	
Suggested A constant*	Interferometry	Ultrasound	Interferometry	Ultrasound
Hoffer Q: pACD	5.35	5.26	5.59	5.35
Holladay 1: Sf	1.60	1.48	1.83	1.57
Barrett: LF	1.78	-	1.86	-
SRK/T: A	118.80	118.59	118.95	118.73
Haigis**: a0; a1; a2	1.36; 0.4; 0.1	1.04; 0.4; 0.1	1.36; 0.4; 0.1	1.13; 0.4; 0.1

* Estimates only: surgeons are recommended to use their own values based upon their personal experience. Refer to our website for updates.

** Not optimized.

INJECTION GUIDELINES


The Medicel Viscoject Bio 1.8 / 2.2 and Accuject 1.8 / 2.0 / 2.1 / 2.2 injection systems are recommended for implanting the FineVision lenses.


These fully single-use systems represent total reliability for safe and effective lens injections. Their compact design with integrated cartridge enables a simple, predictable loading and positioning of the lens.


FineVision Micro F:

Medicel Viscoject Bio 1.8 / 2.2 & Accuject 1.8 / 2.0 / 2.1 / 2.2

Guidelines with Viscoject:

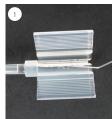

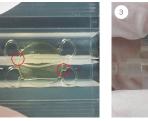

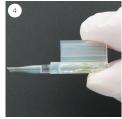






- Apply viscoelastic into the tip and the loading chamber of the injector cartridge.
- 2. Remove the lens from the lens holder. Position the lens into the cartridge in such way that the two haptics with the notches are pointing at 1 and 7 o'clock.
- 3. Exert slight pressure onto the lens optic and make sure that all haptics are inside before further closing the cartridge.
- Close the cartridge and check the position of the lens.
 Once the "click-lock" mechanism engages, the lens is securely loaded.
- 5. Engage the cartridge in the injector.
- 6. Press the injector plunger forward and push the lens into the conical tip of the cartridge. Pull the plunger back a few millimeters and then inject the lens in one continuous motion. For gently implantation, it is not necessary to push the plunger until the end of the cartridge.

FineVision Pod F:

Accuject 2.0 for lens diopters < 25D & Accuject 2.1 / 2.2 for lens diopters \geq 25D Guidelines with Accuject:





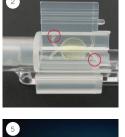





- 1. Apply viscoelastic into the tip and the loading chamber of the injector cartridge.
- 2. Remove the lens from the lens holder. Position the lens into the cartridge in such way that the two haptics with the holes are pointing at 1 and 7 o'clock.
- Exert slight pressure onto the lens optic and make sure that all haptics are inside before further closing the cartridge. Close the cartridge and check the position of the lens.
- 4. Once the "click-lock" mechanism engages, the lens is securely loaded and ready for injection.
- 5. Press the injector plunger forward and push the lens into the conical tip of the cartridge.
- 6. Pull the plunger back a few millimeters and then inject the lens in one continuous motion. For gently implantation, it is not necessary to push the plunger until the end of the cartridge.

Distributed by

FINEVISION_19H1

f in 🖸

PhysIOL sa/nv - Liège Science Park - Allée des Noisetiers 4 - 4031 Liège - Belgium t. +32 (0)4 361 05 49 - f. +32 (0)4 361 05 30 - info@physiol.be - **www.physiol.eu**